




implance.com

AGGRESSOR

IMMEDIATE LOADING

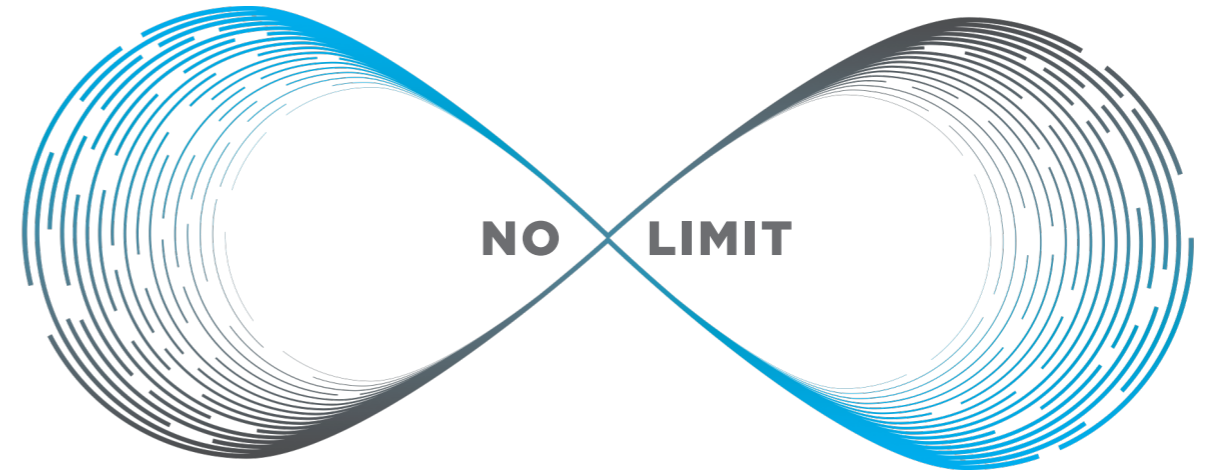

AGS Medikal Ürünleri İth. İhr. Tic. A.Ş.
Katip Salih Sok. No: 111 34718 Koşuyolu / İstanbul / Turkey
T +90 216 546 00 78 F +90 216 546 00 79

implance.com


IMPLANCE
Dental Implant System

AGGRESSOR

Aggressor implant provides unlimited application possibilities.

Aggressor implant offers immediate loading. It does not cause stress in the neck area with its improved aggressive thread design and sharper threads towards the apex.

Aggressor implant body design with increased taper towards its apical, provides high primary stability in the extraction socket and D3-D4 bone type.

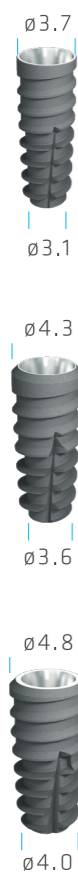

Ø3.7

Length (mm)	Code
10.0	BFAR3710
12.0	BFAR3712
14.0	BFAR3714


Ø4.3

Length (mm)	Code
10.0	BFAR4310
12.0	BFAR4312
14.0	BFAR4314

Ø4.8

Length (mm)	Code
10.0	BFAR4810
12.0	BFAR4812
14.0	BFAR4814

- Aggressive thread design
- Natural root-form
- Sharper threads at the apex
- More contact area with the bone
- Strong primary stability in every bone type
- Immediate implantation and immediate loading protocol

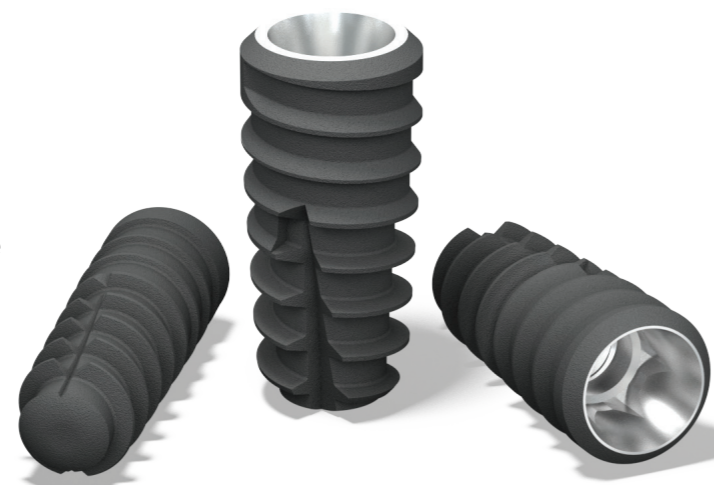


Self-stopper Drills

Self-stopper drills provide great convenience to the users and cancel out the possibility of committing error.

Aggressor implant drills operate at 200-300 rpm, thus preventing necrosis caused by excessive heat. Autogenous bone graft can be collected due to low speed drilling and used for the treatment of the patient.

High primary stability can be achieved with the conical drills designed as per the conical structure of the Implants Aggressor implants.

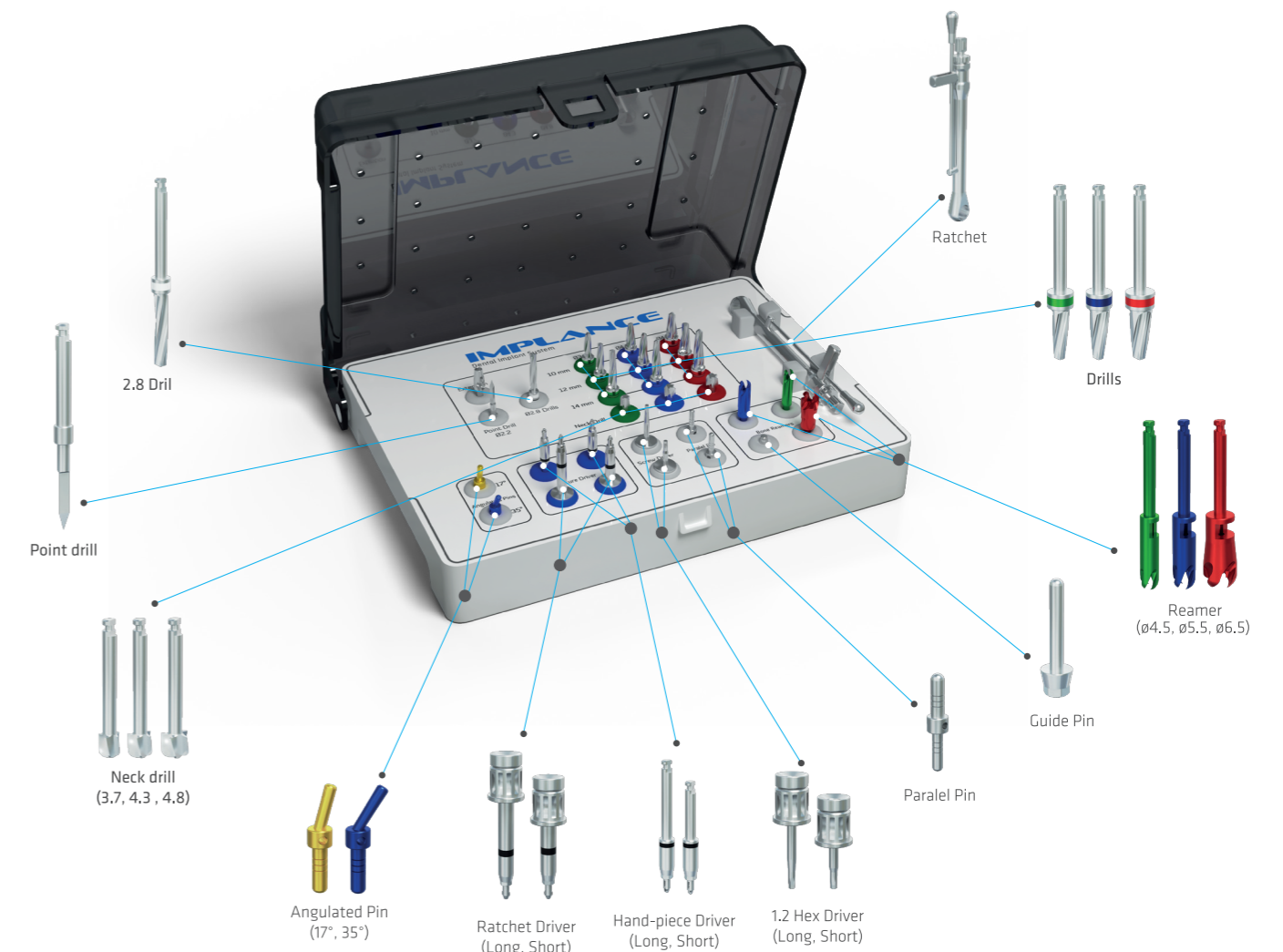

Auxiliary Drills

The ideal exit platform can be created for the abutments by scraping the excess bone in the neck of the implant with the use of Bone Reamer available in three diameters.